AquaFloss Compact

EN – INSTRUCTIONS FOR USE

Product description:

1. Head release button
2. On/Off button
3. "Change mode" button
4. Mode indicator lights
5. Handle
6. Charging station
7. Head
8. Identification ring
9. Suction hose
10. Filter
11. Water reservoir threaded fitting
12. Water reservoir
13. Water reservoir cap

Note: Residual moisture may be found in the unit. This is the distilled water used during product testing.


Instructions:

1. Charging

Plug the charging station into the socket. Place the TrueLife AquaFloss Compact dental water flosser onto the cordless charging stand. The indicator light will flash "CH" while the appliance is charging. When the appliance is fully charged, the light stays on.

Note: If the indicator light flashes "CH" while in use it means the appliance needs to be charged. Charge time is cca 24 hours.

2. Using the AquaFloss Compact dental water flosser

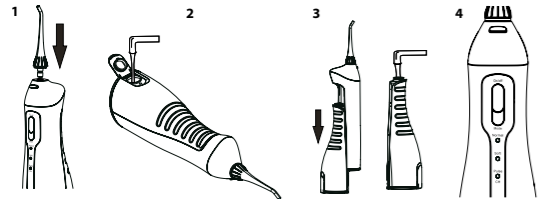- 1 – Attach head to water flosser. It will click into place when correctly inserted.
- 2 – Open water cap and fill reservoir with water no hotter than 40 °C, then close cap and make sure the water is not overflowing.
- 3 – There are two ways to fill the reservoir: the method shown in Fig. 2 and that shown in Fig. 3, where you slide the reservoir down to detach it from the main unit and fill it via the threaded fitting at the top.
- 4 – Press Mode button and choose desired flossing mode. You can change modes while flossing. Modes:
Normal – for regular flossing.
Soft – for gentle flossing of sensitive gums.
Pulse – for flossing and massaging gums.
- 5 – Leave your mouth slightly open to allow water to flow freely from your mouth into the sink. Choose the mode and press the On/Off button to begin flossing. Direct the water jet onto your teeth and gums at an angle of 90°. Slowly move the jet of water around your teeth. Direct water flow onto gums and between teeth. The AquaFloss Compact is suitable for cleaning braces - aim the water jet between the braces brackets.
- 6 – When you want to change heads, press the On/Off button and then the head release button. Remove the head and insert another, if needed.
- 7 – After flossing, turn the appliance off using the On/Off button. Open the water reservoir cap and pour out the remaining water. Turn the appliance on again using the On/Off button to empty the water from the pump. Press the On/Off button to switch off the appliance and close the water reservoir cap. Wipe the appliance with a dry cloth.

Note: To prevent the spread of bacteria, make sure no water remains in the appliance. If you are not intending to use the appliance for some time, make sure it is completely dry before storing.

Package contents:

- 1x TrueLife AquaFloss Compact
- 1x Charging station
- 2x Standard head
- 1x User's manual

Maintenance:

Use only water and pH neutral products when cleaning. Never use corrosive or abrasive cleaners, as they can damage the appliance. Do not use water hotter than 50 °C. Keep appliance away from high temperatures and direct sunlight.

1. Unit

- Wipe away dirt with a towel or cloth and liquid soap and water, then wipe dry with a clean cloth.
- Do not submerge appliance in water.
- Dry after rinsing.

2. Charging station

- Check plug and cable regularly. Clean with cloth.
- Open the charging base cover, then clean and dry the wireless charger. Repeat regularly.

3. Head

- Rinse head with water and wipe with cloth.
- Change heads regularly: at least once every three months depending on wear and tear.

4. Suction hose

- Clean hose with water and dry with a cloth.
- Do not bend, stretch or twist the suction hose.

5. Water reservoir

- Clean reservoir with water and dry with cloth.
- Dry reservoir thoroughly if not intending to use it for some time.

Troubleshooting:

Problem	Possible cause	Solution
TrueLife AquaFloss Compact water flosser not working	Product has just been purchased or has not been used for more than 3 months	Fully charge appliance over 24 hours
Appliance does not work or works for only a short time after full recharge (24h)	Battery nearing end of life cycle	Contact service centre for further action
Water pressure is low	Head is not clean / Filter is clogged with dirt / Head is misshapen	Insert new head / Clean filter
No water coming out of head	Reservoir is empty / Unit is tilted too much	Fill reservoir / Put unit in correct position (vertical)

Manufacturer's Liability:

TrueLife products are subject to warranty for manufacturing defects. The obligations of the company, elem6 s.r.o., in the case of such defects is to repair or replace the defective component or product under warranty for an equivalent component or product, at the manufacturer's discretion. Transfer of ownership, alteration or unauthorised manipulation of the product exceeds the scope of this warranty, even in cases where defects, negligence and other manufacturing errors result in loss or damage. The warranty does not extend to damage due to use, accident or normal wear and tear. The company, elem6 s.r.o., is not liable for any accidents, injuries, death, loss or other claims related to or resulting from the use of this product. Under no circumstances does the company, elem6 s.r.o., accept liability for incidental or consequential damage related to or resulting from the use of this product or any of its components. Return or exchange of spare parts and/or products may incur fees for transportation, handling, replacement and/or refilling.

This product should be stored out of the reach of children. Children should use it only under adult supervision.

SAFETY WHEN USING THE APPLIANCE:

1. Carefully read the enclosed safety instructions.
2. For correct handling of the product, refer to the latest version of the manual, which is available for download at www.truelife.eu
3. Elem6 s.r.o. will not be held responsible for misprints and clerical errors.

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